

Medela's Invia[®] NPWT Portfolio

FACT SHEET

Therapy Confidence, Easy Application, Simplified Management

- **Dynamic** pressure management that provides reliable negative pressure delivery from pump to wound site¹
- **Proven** pump performance and versatile portfolio³
- **Easy** dressing use and application²
- **Secure** canister and tubing connection²




Medela simplifies the delivery and management of NPWT making it easier for healthcare professionals to help their patients heal.

Invia NPWT systems actively control set pressure at the wound site and adapt the flushing cycle to the amount of fluid in the tubing for optimal system performance. When using the Medela dressing kits with FitPad in conjunction with the new Invia Negative Pressure Wound Therapy (NPWT) Systems, healthcare professionals can be confident that negative pressure is delivered¹.

Invia® Liberty™ NPWT System

- Invia Liberty offers active exudate management and pressure control from the pump to the wound site
- The pressure options reach from –40 to –200 mmHg including standard –125 mmHg
- Invia Liberty has an on-screen air leak indicator and new custom air leak modes for normal and large wounds
- Easy and secure tubing attachment is possible with the new Quick-connector
- The canisters are easily attached and detached with a one-touch canister release button
- Healthcare professionals have the choice of 300 ml and 800 ml canisters with solidifier

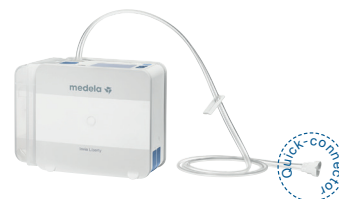

Invia® Motion™ NPWT Systems

Invia Motion delivers the convenience of a personal pump while offering clinical flexibility with a choice of pressure settings and therapy modes.

- Providing NPWT in the hospital, in transition, and at home
- Invia Motion offers a clear user interface supported by 4 buttons and a digital display
- It comes with a discrete carrying case, battery autonomy of 10 hours and compact design for comfort and mobility
- The preset pressure settings reach from –40 to –175 mmHg (–125 mmHg default setting)
- Easy and secure tubing attachment is possible with the new Quick-connector
- With Invia Motion clinicians and patients have the choice of therapy units to meet clinical and mobility needs



Invia® NPWT Dressings

- The Invia NPWT Dressings support active exudate management and pressure control from the pump to the wound site
- The Dressings allow an easy application of transparent film and external suction interface (FitPad)
- The low profile design and flexible wings of the suction interface (Invia FitPad) allow appropriate placement in difficult anatomic locations
- Easy and secure tubing attachment is possible with the new Quick-connector
- Medela offers a streamlined foam and gauze dressings portfolio



Invia Silverlon NPWT Dressing

Silver ions delivered in the dressing when activated by moisture provide a barrier to bacterial penetration which can help reduce infection.

- Can be cut or overlapped to fit any wound size.
- Non-adherent for comfortable dressing changes¹.



For specific NPWT indications, contraindications, warnings and cautions, consult the appropriate Invia NPWT instructions for use.

Obtain more information about Medela's NPWT portfolio at www.medela-healthcare.us info-healthcare@medela.com or 1-877-735-1626.

¹ Performance data on file at Medela AG

² Usability Validation data on file at Medela AG

³ Invia Liberty system is available in the US market since 2009 and over 100'000 patients have been treated with Invia Liberty (in the US) and Avance devices (in Europe)